Gebruiksaanleitung Manual Mode d'emploi

D

Seite 2

GB

Page 4

F

Page 6

Alpha[®] Drehscheibe
Alpha[®] Turning Disc
Alpha[®] Disque Pivotant

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1. Specifications and Scope of Supply

Item no.: PM-5080
Dimensions: Ø 40 cm
Material: ABS coated plastic fabric, latex free
Colour: black
ILN: 4250355350800
Storage: -5 – 40° C
 protect from sunlight


Scope of supply:

- Turning disc
- Manual

Made in Germany



max.
200 kg

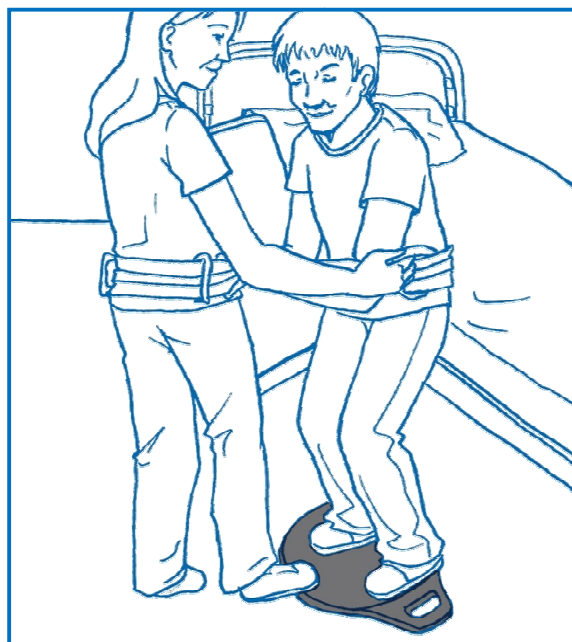2. Use / Placement

The turning disc is used to transfer a patient capable of sitting who can take on some weight with their feet. It is also ideal for use in conjunction with a sliding mat, raising aid, sliding board or belt.

This product is suitable for re-use.

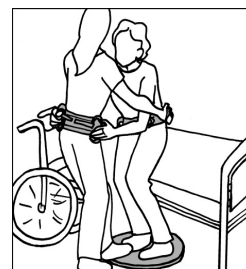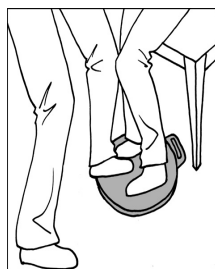3. Maintenance / Care

The product does not require any special maintenance. You can use a normal household cleaner and wipe it down with a moist cloth regularly for reasons of hygiene. A wipe-down disinfection can be performed with the usual disinfectant solutions and isopropyl alcohol. There is also some care information on the label.

4. Warranty

- Petermann gives a warranty that the product is devoid of material defects and manufacturing defects.

- The warranty is valid for a period of **3 years**. It does not cover any defects which are due to improper use.
- Petermann cannot be held liable for any accidental or indirect injury or damage.

5. Visual and Function Inspection / Safety Information / Re-use

Repositioning products are also subject to natural aging. Check the equipment regularly and carefully to make sure it is working properly and especially before re-use. Make sure that the product is in good condition and that the material is not damaged. Make sure that the material resists the load and does not break. The special material of the turning disc allows for easy turning. The best way to test these properties is by comparing the product with a new product. For re-using the product it must be disinfected before, see 3.

WARNING:

- ***Faulty products must no longer be used.***
- ***If the patient's mobility is restricted, it is essential to supervise the motions of transfer by providing trusted help.***
- ***Make sure the patient's feet are properly and entirely placed on the turning disc and that the turning disc is lying firmly on the floor.***
- ***Make sure the turning disc does not constitute an obstacle in the room and no one can trip over it.***
- ***Carefully check the transfer surfaces, making sure they are stable and brakes are applied (e. g. wheelchair, bed).***

6. Product Labeling

The Petermann Alpha[®] turning disc is a medical product, class I and complies with the MDD directive 93/42/EEC.

Petermann GmbH is certified in accordance with the directive DIN EN ISO 13485 by the German TUV.

7. Disposal Information

For the turning disc there are no special disposal rules; it can be disposed of in accordance with national and local regulations as necessary.

If you have any questions about the use of our products, please call our service number:

+49 (0) 9868 9339-0

Petermann GmbH • Ziegelhaus 12 • D-91601 Dombühl
 Postal address: Petermann GmbH • Waldweg 3 • D-91601 Dombühl
 Phone: +49 (0) 9868 9339-0 • Fax: +49 (0) 9868 9339-20
 E-mail: petermann@pm-med.de • www.pm-med.de

"Information subject to technical modifications to our products and deviations in colour or printing errors"

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Petermann GmbH est certifié par le TUV conforme à la directive DIN EN ISO 13485.
Adresse postale: Petermann GmbH • Waldweg 3 • D-91601 Dombühl